

# Esophageal Dilatation Systems



# M-Flex® Blue Silicone Bougies



## There's nothing more flexible than an M-Flex.®

Our silicone bougies feature added flexibility that some users prefer. Tungsten-filled, there's nothing quite like the M-Flex.®

### Benefits

- Ultra-flexible tapered tip (Maloney)
- Completely radiopaque from distal tip to proximal end
- American and European depth markings
- Mercury-free
- **Autoclavable for quick reprocessing**

### Features

- 21 sizes: 20 FR to 60 FR (in 2 FR increments)
- Choose a full set, mini set or individual pieces
- Easy disposal
- No charge trade-in program for proper disposal of expired mercury-filled bougies (with purchase of new product)




Store and transport M-Flex® Bougies in our storage case which holds up to two full sets.

### M-Flex® Dilatation Systems

Set	Part No.	Contents
Maloney (tapered) 21-piece full set	1206-00	21 bougies, 20–60 FR
Maloney (tapered) 10-piece mini set	1206-10	10 bougies, 36–54 FR
Hurst (blunt) 21-piece full set	1207-00	21 bougies, 20–60 FR
Hurst (blunt) 10-piece mini set	1207-10	10 bougies, 36–54 FR
Bougie Storage Case	1210-01	Heavy-duty storage for up to 46 bougies. 32"L x 7"W x 4"H



Our storage case allows easy identification of bougie size. Removing plexiglass window provides easy access.

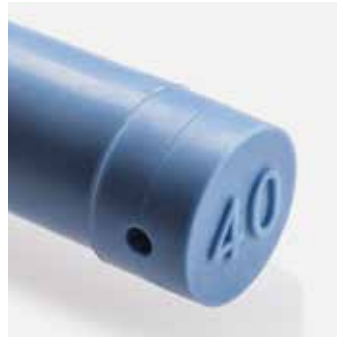

End caps display French size.



Extra flexible tip.



Top: Maloney tip. Bottom: Hurst tip.

### M-Flex® Blue Silicone Bougies

FR (French) Size	(mm)	Maloney (tapered) Tip Part No.	Hurst (blunt) Tip Part No.
20	6.6	1206-20	1207-20
22	7.4	1206-22	1207-22
24	8	1206-24	1207-24
26	8.6	1206-26	1207-26
28	9.3	1206-28	1207-28
30	10	1206-30	1207-30
32	10.7	1206-32	1207-32
34	11.2	1206-34	1207-34
36	12	1206-36	1207-36
38	12.5	1206-38	1207-38
40	13	1206-40	1207-40
42	14	1206-42	1207-42
44	14.5	1206-44	1207-44
46	15.2	1206-46	1207-46
48	16	1206-48	1207-48
50	16.5	1206-50	1207-50
52	17.2	1206-52	1207-52
54	18	1206-54	1207-54
56	18.5	1206-56	1207-56
58	19.2	1206-58	1207-58
60	20	1206-60	1207-60

Hurst Bougie

Maloney Bougie

American Depth Markings

European Depth Markings



# WeightRight® Mercury-Free Bougies



## Heavy-weight control and convenience combined.

Our tungsten-filled bougie mimics the weight and feel of mercury-filled bougies without the concerns of mercury usage and disposal. It's the weight you expect, and the convenience you demand.

### Benefits

- Completely radiopaque from distal tip to proximal end
- Flexible material does not attract particulate
- American and European depth markings
- Mercury-free

### Features

- 23 sizes: 16 FR to 60 FR (in 2 FR increments)
- Choose a full set, mini set or individual pieces
- Easy disposal
- No charge trade-in program for proper disposal of expired mercury-filled bougies (with purchase of new product)



Store and transport WeightRight® Bougies in our storage case.

### WeightRight® Dilatation Systems

Set	Part No.	Contents
Maloney (tapered) 23-piece full set	1208-00	23 bougies, 16–60 FR
Maloney (tapered) 10-piece mini set	1208-10	10 bougies, 36–54 FR
Hurst (blunt) 23-piece full set	1209-00	23 bougies, 16–60 FR
Hurst (blunt) 10-piece mini set	1209-10	10 bougies, 36–54 FR
Bougie Storage Case	1210-01	Heavy-duty storage for up to 46 bougies. 32"L x 7"W x 4"H





Our storage case allows easy identification of bougie size. Removing plexiglass window provides easy access.



Hang bougie vertically using the hole in the end cap.



Top: Maloney tip. Bottom: Hurst tip.

Hurst Bougie

Maloney Bougie

American Depth Markings

European Depth Markings



### WeightRight® Mercury-free Bougies

FR (French) Size	(mm)	Maloney (tapered) Tip Part No.	Hurst (blunt) Tip Part No.
16	5.3	1208-16	1209-16
18	6	1208-18	1209-18
20	6.6	1208-20	1209-20
22	7.4	1208-22	1209-22
24	8	1208-24	1209-24
26	8.6	1208-26	1209-26
28	9.3	1208-28	1209-28
30	10	1208-30	1209-30
32	10.7	1208-32	1209-32
34	11.2	1208-34	1209-34
36	12	1208-36	1209-36
38	12.5	1208-38	1209-38
40	13	1208-40	1209-40
42	14	1208-42	1209-42
44	14.5	1208-44	1209-44
46	15.2	1208-46	1209-46
48	16	1208-48	1209-48
50	16.5	1208-50	1209-50
52	17.2	1208-52	1209-52
54	18	1208-54	1209-54
56	18.5	1208-56	1209-56
58	19.2	1208-58	1209-58
60	20	1208-60	1209-60

# SafeGuide® Over the Guidewire Esophageal Dilators



## Flexibility and control for exacting procedures.

At times, there is need for greater accuracy. That's why there's SafeGuide®, a remarkably flexible over the guidewire solution that ensures exact placement of the dilator. Position the guidewire with the assistance of an endoscope and direct the SafeGuide® dilator where you need it.

### Benefits

- Completely radiopaque from distal tip to proximal end
- Flexible material does not attract particulate
- American and European depth markings
- Dilators are reusable

### Features

- 16 sizes: 15 FR to 60 FR (in 3 FR increments)
- Choose a full set, mini set or individual pieces
- Easy to clean with optional brushes and adapters
- With proper care, reusable SafeGuide® Dilators provide years of service



Store and transport SafeGuide® Dilators in a convenient storage case. A labeled slot holds each dilator secure. The case holds a full set of 16 dilators.

### SafeGuide® Dilatation Systems

Set	Part No.	Contents
16-piece full set	1214-00	16 dilators, 15–60 FR, 1 storage case, 1 marked spring tip guidewire, 1 reusable cleaning brush, 25 disposable cleaning adapters
7-piece mini set	1214-10	7 dilators (15, 21, 27, 33, 39, 42, 45 FR), 1 storage case, 1 marked spring tip guidewire, 1 reusable cleaning brush, 25 disposable cleaning adapters
Custom mini set	Contact your sales representative (Available in US only)	7 dilators of your choice, 1 storage case, 1 marked spring tip guidewire, 1 reusable cleaning brush, 25 disposable cleaning adapters



The single-piece construction of the SafeGuide® Spring Tip Guidewire resists breakage and flexes easily.

## SafeGuide® Over the Guidewire Esophageal Dilators

FR (French) Size	(mm)	Part No.
15	5	1214-15
18	6	1214-18
21	7	1214-21
24	8	1214-24
27	9	1214-27
30	10	1214-30
33	11	1214-33
36	12	1214-36
39	13	1214-39
42	14	1214-42
45	15	1214-45
48	16	1214-48
51	17	1214-51
54	18	1214-54
57	19	1214-57
60	20	1214-60

## Accessories



Storage case indicates size.



Clearly marked reusable and single-use cleaning brushes.

## SafeGuide® Accessories

Set	Part No.	Description	Qty/Case
SafeGuide® Storage Case	1214-01	Holds up to 16 SafeGuide® Dilators. 36"L x 10"W x 4"H	1
Reusable Marked Spring Tip Guidewires	1214-02	Guides the dilator. 210cm	2
Disposable Marked Spring Tip Guidewires	1214-02D	Guides the dilator. 210cm	6
Disposable Cleaning Brushes	1214-03D	Cleans the inner lumen of the dilator	3
Reusable Cleaning Brushes	1214-03R	Cleans the inner lumen of the dilator	3
Disposable Cleaning Adapters	1214-04	Used to flush the inner lumen of the dilator	25



Insert cleaning brush into the distal end and push through to the other side for complete inner lumen cleaning.



Disposable cleaning adapter.

SafeGuide® Dilator



## Is it time to replace your esophageal bougies or wire-guided dilators?

### Replace bougies when you notice:

- End cap separating from bougie
- Cuts, cracks or nicks
- Flaking or peeling of the outer coating
- Faded or missing markings
- Expired or missing expiration dates

### Replace wire-guided dilators when you notice:

- Loss of flexibility
- Cuts, cracks or nicks
- Stretching in the inner lumen
- Faded or missing markings
- Curved dilator or curled tapered tip
- Expired

## Is it time to replace your esophageal dilator guidewire?

### Replace guidewires when you notice:

- Kinked or bent wire
- Kinked, bent, or stretched out spring tip

## Is it time to replace your esophageal dilator cleaning brushes?

### Replace brushes when you notice:

- Kinked or bent wire
- Brush bristles falling out

**Diversatek™**  
Healthcare

[DiversatekHealthcare.com](http://DiversatekHealthcare.com)

P 800.558.6408 or 414.265.7620

#### Corporate Headquarters

102 East Keefe Avenue  
Milwaukee, WI 53212 USA

orders@diversatekhc.com  
sales@diversatekhc.com

#### Technical Research & Training Center

9150 Commerce Center Circle, Suite 500  
Highlands Ranch, CO 80129 USA

Learn more from  
your personal product  
representative at  
**800-558-6408**