

Giving you more choices with Maquet Meera CL

Just what you need – no more, no less

Every hospital needs solid, functional OR tables, but not every OR needs a table with every feature. That's why the Maquet Meera family now includes Maquet Meera CL, ideal for versatile use or general surgery.



Explore the benefits of Maquet Meera CL

Excellent ratio of price to performance

Not all hospitals require motorized leg adjustments. Maquet Meera CL is built on the same components as the rest of the Maquet Meera line, but on a simplified table base with a mechanical locking system and manual leg movement. This gives you the most important of Maquet Meera features, but at a reduced price.

Ergonomic handling

The new Smart Control with graphic touchscreen makes setup easy. It continuously receives and displays real-time table configuration information, which saves time, enhances safety, and simplifies documentation.

Stability for the largest patients

The robust base ensures stability and safety for patients weighing up to 454 kg / 1,000 lbs. For the complex demands of bariatric surgery, Maquet Meera CL offers an unrestricted capacity of 250 kg / 551 lbs for any table setup.

Wide range of compatible accessories

All Maquet Meera Tables feature the same interfaces as past Getinge Mobile OR Tables to ensure compatibility with existing Getinge accessories. This increases the flexibility of existing equipment to extend the life of your existing OR table investment.

Longitudinal shift improves access for better outcomes


The motorized longitudinal shift delivers increased access for the surgeon, while simultaneously improving accessibility for the C-arm, improving intraoperative diagnostics and treatment. 360° imaging is possible when paired with carbon fiber accessories.

Maneuverability

The special design of the castors enables intermediate operation, straight-line travel, and free-wheeling. The ability for lateral movement is ideal for outpatient surgery centers.

Construction features

- Stable base construction with four double swivel castors for easy movement and maneuvering (base can be locked via foot pedal)
- Cover for the override panel made of GFR composite plastic, resistant to impact, breakage and disinfectants
- Column casing made of chrome-nickel steel
- Identical interfaces on normal and on reverse side

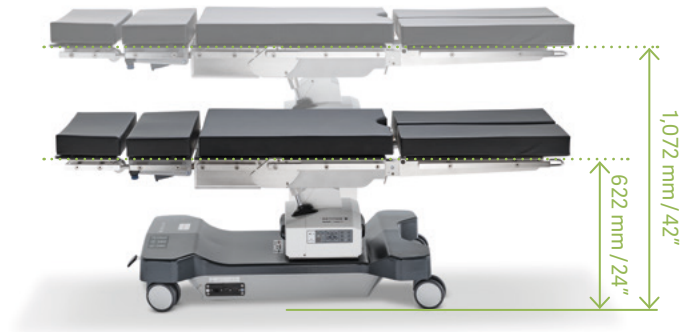

Electrical specifications

- Specially designed rechargeable batteries, with a single charge lasting about one week in the operating room
- Electronic charge monitoring, with visual and aural indicators
- Batteries recharged from the mains supply
- 100–240 V AC (adjustable), 50–60 Hz, over mains cable

Smart Control features

- Display of table set up information
- Memory function for up to 10 table top positions
- Charging status of column and hand control
- Service and error messages
- Pre-programmed positions: flex, reflex and back-horizontal
- Table lock for blocking movements
- Backlight function

Characteristics Maquet Meera CL 7000.01X0

Weight	218 kg / 481 lbs
Length of OR table without accessories	765 mm / 30"
Maximum overall load (patient plus accessories)	454 kg / 1,000 lbs in Normal, 250 kg / 551 lbs in Reverse
Trendelenburg / Reverse Trendelenburg	25° / 35°
Lateral tilt	20°
Lower back plate	+75° / -40°
Leg plate (manual)	-90°
Motorized longitudinal shift	310 mm / 12"

This document is intended to provide information to an international audience outside of the US.
Maquet GmbH · Kehler Str. 31 · 76437 Rastatt · Germany · +49 7222 932-0 · info.sales@getinge.com