

Ninjo™ Flusher Disinfecter

Attractive modern design, easy operation and excellent disinfection performance




Choose from Ninjo™ colourful polymeric or stainless steel finish, and either free standing, wall-mounted or under-counter models, to suit your sluice room needs and preferences.

- **Digital A₀ value** and remaining time display
- **Built-in cooling fan** dries containers quickly and cools them with HEPA-filtered (H13) air.
- **Reduce risk of cross-contamination**
- **Two programmes available:** economy and normal, with an additional three available upon request.
- **Flexible design options**

Ninjo™ cleans receptacles with a series of high pressure water cycles delivered through strategically located fixed and rotating nozzles, followed by a thermal steam disinfection phase.

Specifications

PRODUCT INFORMATION		
Ninjo	Free Standing	Under Counter
Weight	90 kg (198 lbs)	70 kg (155 lbs)
Width	450 mm (17 ¹¹ / ₁₆ "	450 mm (17 ¹¹ / ₁₆ "
Depth	580 mm (22 ³ / ₄ "	580 mm (22 ³ / ₄ "
Height	1320 mm (51 ⁷ / ₈ "	870 mm (34 ³ / ₁₆ "
Max. external temperature	40° C (100° F)	
Noise level	53 dB (A) accordingly to EN ISO 3747:2000	
Electrical connection	Contact your Arjo representative	
Power Consumption	0,16 kWh/cycle	
Water Consumption	Economy 11 liters/cycle ±10 % Normal 18 liters/cycle ±10 %	
Water Connection	Pressure 1-8 Bar (14-116 PSI) 15 mm (1/2"), Rate of flow >0,3 l/s	
Drain Connection	Wall ø 90 or 110 mm (3 1/2 or 4 1/4") Floor ø 90 or 110 mm (3 1/2 or 4 1/4")	
Approvals	MDD 93/42/EEC, Third party certificate for EN ISO 15883-1 and 15883-3	
Other	Dimensions and specifications of the Ninjo matching cabinets available on request	

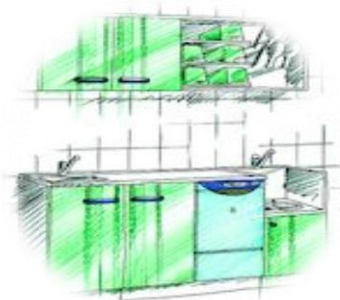

At Arjo, we are committed to improving the everyday lives of people affected by reduced mobility and age-related health challenges. With products and solutions that ensure ergonomic patient handling, personal hygiene, disinfection, diagnostics, and the effective prevention of pressure injury and venous thromboembolism, we help professionals across care environments to continually raise the standard of safe and dignified care. Everything we do, we do with people in mind.