

Customer:

Reference:

INSPECTION LAMP, STANDARD

PRODUCT DESCRIPTION


Inspection lamp with integrated magnifier for professional use, especially suitable wherever demanding work is required. With its effective, anti-dazzle light, without annoying heat, it is the ideal product for checking delicate instruments. The position of the inspection lamp can be easily adjusted with one hand, enabling staff to find the optimal height and angle for the magnifier. The high quality lens magnifies x1.75 and the size of the viewing field is 127mm. The circular shape of the 22W fluorescent lamp ensures that the 1350 lumen light is distributed evenly. A dust cover is mounted on the magnifier lens to protect it from dust and from inadvertently acting as a burning glass. The arm of the inspection lamp is 980 mm and is supplied with a clamp for fixation to the work surface.

APPLICATION

Designed for work with visual inspection at the CSSD or other smaller sterile units in hospitals and clinics.

KEY FEATURES

- Ergonomic
- Practical
- Durable

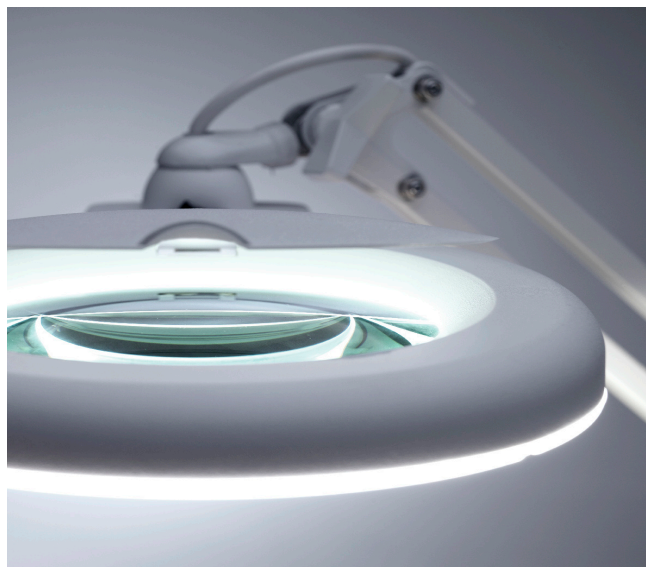

STANDARDS & CODES

Each lamp is built in accordance with the following EU safety norms:

- CE-marked
- Compliant with RoHS directive
- Compliant with WEEE directive

PRINCIPLE OF OPERATION

Operation is extremely simple. The lamp is made ready for use by adjusting it to any desired position and switching it on.

TESTING

All inspections lamps have been tested in accordance with Getinge demands.

TECHNICAL DATA

Mark your selection.

Art. No.	Description	Voltage
561 95 26-01 <input type="checkbox"/>	Inspection lamp with EUR plug	220-240V / 50Hz
561 95 26-02 <input type="checkbox"/>	Inspection lamp with UK 3-pin plug	220-240V / 50Hz