



DetachaTip® III System Multi-use Endosurgery Instruments

Only the DetachaTip® III System offers this combination of performance, reliability and cost efficiency.



3DHD
VISION SYSTEM



UNIVERSAL PLUS®
SUCTION / IRRIGATION



ENTRIPORT*
TROCAR SYSTEM



To learn more about these and other innovative products, call **800-448-6506** or visit **ConMed.com**.

DETACHATIP® III SYSTEM

NEW



LONG FENESTRATED

NEW



CLAW

NEW



BLUNT NOSE

Internal Ratchet Design

- Convenient location for standard use or palming; right or left handed operation.
- Audible and tactile feedback.
- Includes innovative intermittent on/off functionality.

Note: When using any scissor shaft, the ratchet is automatically disabled.

Unmatched performance and cost efficiency

The DetachaTip® III System provides complimentary balance in a reusable handheld instrument system. It combines precise, reliable, and comfortable performance with dramatically reduced procedural costs. Only DetachaTip III offers the optimal blend of performance and cost efficiency.

Ergonomic comfort

The DetachaTip III handle is ergonomically redesigned to accommodate a variety of surgeon hand sizes, is textured for a better grip and weight balanced for greater stability with positive feedback while minimizing hand fatigue. The number of reprocessing cycles (sterilization cycles) has been significantly increased which extends product life resulting in significant cost efficiencies.

Electrosurgical Monopolar Cautery Post


Compatible with monopolar electrosurgical cables.

Shaft Lock/Button Release

- Convenient location to prevent accidental shaft release.
- Provides convenient right or left handed operation to release shaft from handle.

Thumb and Finger Loops

- Comfortable fit for many hand sizes.
- Patented designs for physiologic positioning reduces hand fatigue during procedure.

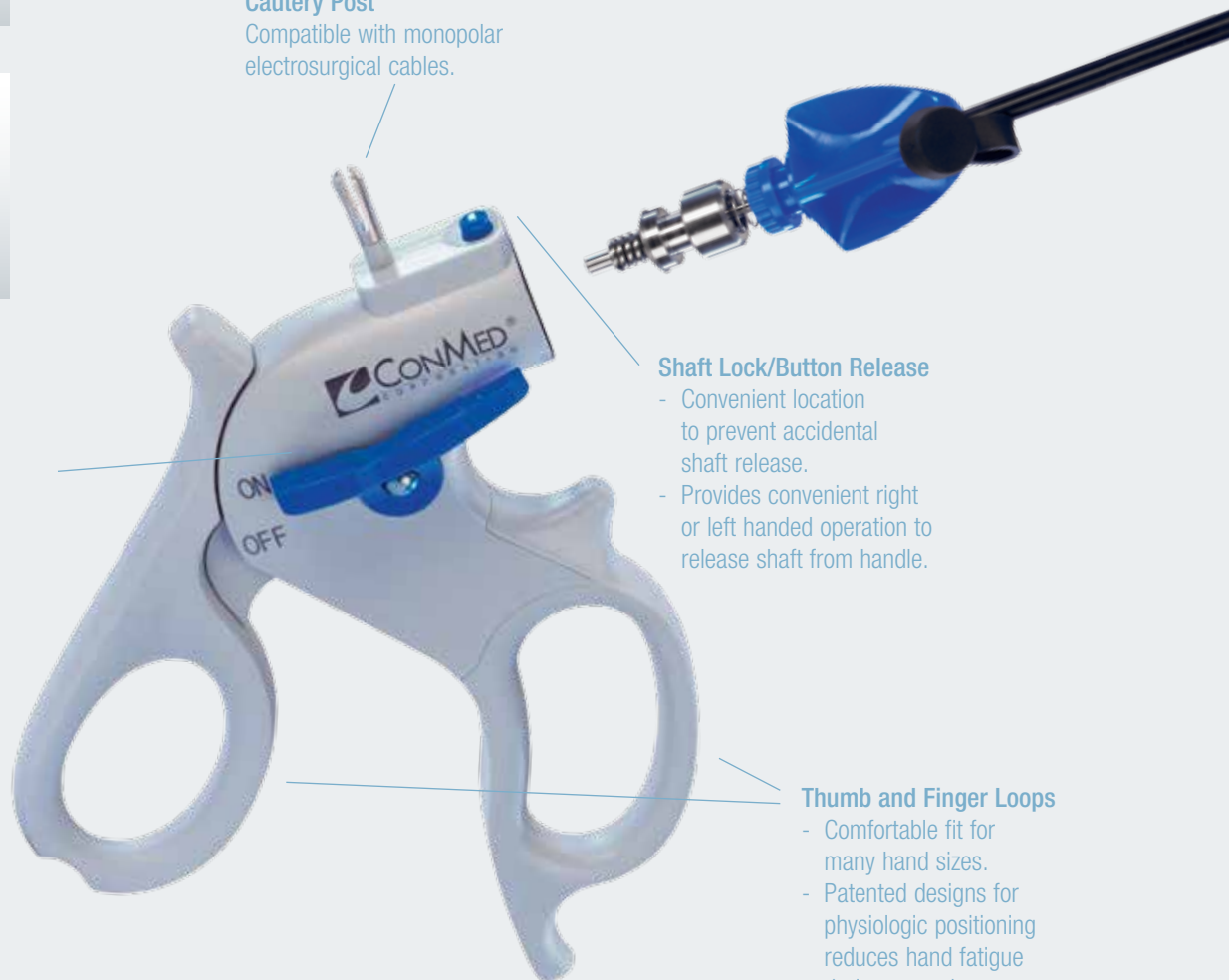

NEW DetachaTip® III Sterilization Trays

DetachaTip® III Instrument Sterilization Trays are designed specifically for the DetachaTip III shafts and handles. The DetachaTip III PLUS Sterilization Tray accommodates the DetachaTip shafts and handles in the top tray, and other ConMed multi-use devices in the bottom of the tray.



Advanced Rotation Mechanism

24 clicks provide very precise control when rotating shaft.

Rotator

Accommodates entire fingertip for comfortable, precise control of tip rotation.

Shaft

New composite shaft material does not conduct electricity, increases durability, minimizes risk of shock.

New Luer Lock Cap and Port

Innovative design allows high flow lumen flushing for manual or automated systems for safe, fast, effective cleaning. No need to “take apart” shafts.

Multi-Use Scissors

S METZENBAUM CURVED

S MINI-METZENBAUM CURVED

Multi-Use Graspers

G STRAIGHT FENESTRATED (DUCKBILL)

G BABCOCK 5MM

G ALLIS

G ENDOWEAVE™

Multi-Use Dissectors

D TAPERED (DOLPHIN)

D MEEKER RIGHT-ANGLE

D MARYLAND CURVED



525 French Road, Utica, NY 13502

DETACHATIP® III
SYSTEM

DetachaTip® III System Ordering Information

Description	Unit	Catalog Number
DetachaTip® III Scissors		
Curved Metzenbaum, 33cm x 5mm	1/ea	3-1003
Curved Metzenbaum, 43cm x 5mm	1/ea	3-4301
Curved Mini-Metzenbaum, 33cm x 5mm	1/ea	3-1004
Curved Mini-Metzenbaum, 43cm x 5mm	1/ea	3-4304
DetachaTip III Graspers		
Straight Fenestrated Grasper, 33cm x 5mm	1/ea	3-1008
Straight Fenestrated Grasper, 43cm x 5mm	1/ea	3-4307
Endoweave™ Grasper, 33cm x 5mm	1/ea	3-1028
Endoweave™ Grasper, 43cm x 5mm	1/ea	3-4328
Allis Grasper, 33cm x 5mm	1/ea	3-1019
Allis Grasper, 43cm x 5mm	1/ea	3-4319
5mm Babcock Grasper, 33cm x 5mm	1/ea	3-1005
5mm Babcock Grasper, 43cm x 5mm	1/ea	3-4305
Long Fenestrated Grasper, 33cm x 5mm	1/ea	3-1006
Long Fenestrated Grasper, 43cm x 5mm	1/ea	3-4306
Claw Grasper, 33cm x 5mm	1/ea	3-1011
Claw Grasper, 43cm x 5mm	1/ea	3-4311
DetachaTip III Dissectors		
Maryland (Curved) Dissector, 33cm x 5mm	1/cs	3-1009
Maryland (Curved) Dissector, 43cm x 5mm	1/cs	3-4308
Meeker (Right Angle) Dissector, 33cm x 5mm	1/cs	3-1017
Meeker (Right Angle) Dissector, 43cm x 5mm	1/cs	3-4317
Tapered (Dolphin) Dissector, 33cm x 5mm	1/cs	3-1018
Tapered (Dolphin) Dissector, 43cm x 5mm	1/cs	3-4318
Blunt Nose Dissector, 33cm x 5mm	1/cs	3-1012
Blunt Nose Dissector, 43cm x 5mm	1/cs	3-4312
DetachaTip III Handle		
Standard Handle with Ratchet	1/cs	3-1010
DetachaTip III Sterilization Trays		
DetachaTip III Tray	1/cs	3-4329
DetachaTip III PLUS Tray	1/cs	3-4343
DetachaTip III Accessory		
Luer Lock Cap	5/box	CD203

Ask about our **NEW 3DHD Laparoscopic Visualization Technology** – bringing depth and clarity to surgery.

Customer Service: 1-800-448-6506
International Sales: 1-315-797-8375
Fax: 1-800-438-3051

ConMed.com
info@mail.conmed.com