

AUTOCLAVE
PRIMA



Autoclave Prima

TECHNICAL EVOLUTION, GREAT DESIGN

Our Galbiati Lab, 100% Italian

Department of “research and development”, created PRIMA, Class B **steam sterilizer** with an uninitiated design sense of aesthetic cleanliness combined with technology state-of-the-art and easy to use.

Solid, easy to use with an uncommon elegance, it's the new born of the Galbiati family, the first Italian production company of dental units.


TECHNICAL SPECIFICATION

- **Rated voltage:** a.c.220V-230V, 50Hz
- **Rated power:** 18L 1500VA, 23L 1700VA
- **Material:** stainless steel (for medical)
- **Max. work pressure:** 2.5 bar Min. work pressure: -0.9 bar
- **Max. temperature:** 145 °C
- **Chamber volume:** 18L (Φ245×320mm) 23L (Φ245×450mm)
- **Loading size:** 18L (198×204×285) 23 L(198×204×385)
- **Max. loading weight:** 18L (3.07kg/cm²)
23L (3.21kg/cm²)
- **Working pressure/temperature:** 1.10~1.30bar /
121°C~122°C; 2.10~2.30bar / 134°C~135°C
- **Water volume for one cycle:** 0.16L (min) 0.18L(max)

PERFORMANCES

- Fast cycles (26min for HOLLOW 134°C)
- Limited water consumption
- Automatic water download
- 100% stainless steel
- Electronic water quality control
- Easy to position system
- GKE helix conformity test
- Integrated printer and usb data saving systems
- 24 months warranty

DIMENSIONS

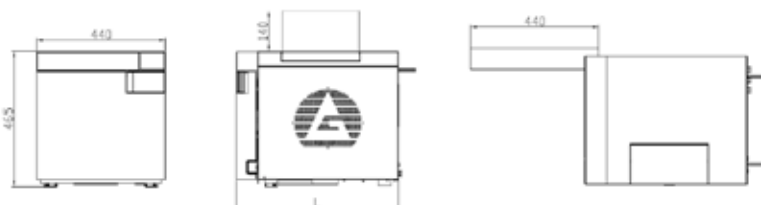

Item	18L Model	23L Model
Size of Closed Door	L615×465×440	VL715×465×440
Size of Open Door	L1005×605×440	L1105×605×440

Button-control panel

Simple and intuitive

Clean watertank


Stainless steel tank, indestructible

Display

Cycle diagram display for live work control

Chamber

4 storey tray chamber, usable also for surgical kits

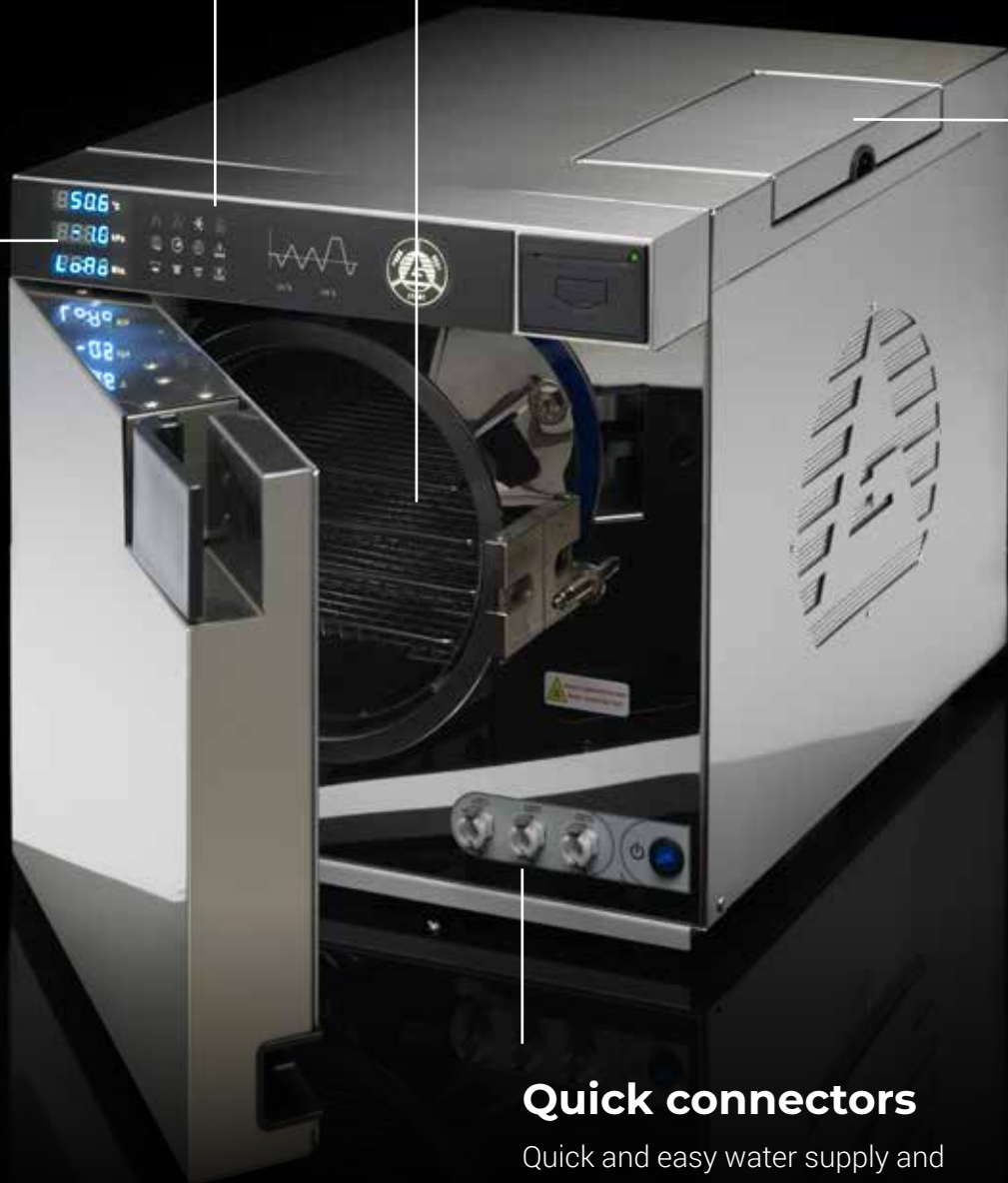

Quick connectors

Quick and easy water supply and drainage connections



Main PCB board

Easy to reach PCB board for maintenance

Galbiati

EXPERIENCE AND QUALITY SINCE 1935

“When a product does not work how it should, or our customers don’t like it, for us is a challenge. We feel that from our factory only good products, functional and accessible to all must be produced.

Our daily commitment is to satisfy our customer at 360° and be present before, during and after the purchase of Galbiati products, with the same passion and dedication”.



www.galbiati.com
info@galbiati.com