

Visible Staple Counter for constant indication of staple supply and for efficient use.

Arrow Guide for accurate and consistent placement.

Open Nose Piece and Low Profile for clear view of staple placement.

REFLEX® TL

REFLEX® ONE

CONMED REFLEX® SKIN STAPLERS

Unequalled Quality, Unequalled Control

- Two stapler styles — Reflex® One and Reflex® TL
- Two stapler sizes — Regular and Wide
- Confirmation of complete staple formation
- Excellent visibility
- Visible staple counter for more efficient use


Ergonomically Designed to help reduce patient discomfort.

REFLEX® EXTRACTOR


CONMED®
ENDOSURGERY

CONMED Reflex® Skin Staplers

Delivering The Performance You Demand!

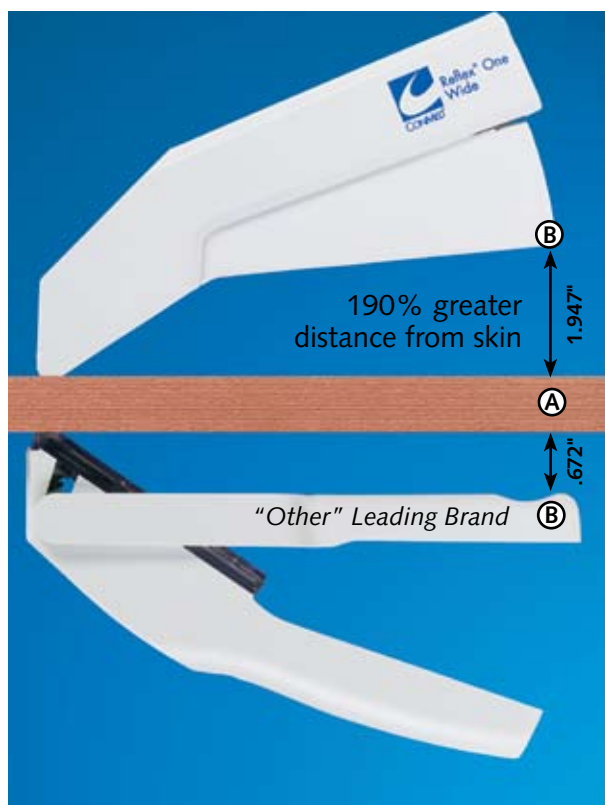

CONMED's Reflex® One and Reflex® TL ergonomic design gives you more area between the handle and the skin (B>A). This keeps the nose of the stapler parallel to the skin allowing proper staple formation, good tissue proximation, and less chance of the onset of edema. When edema occurs, tissue grows over the staple, making for difficult and painful staple extraction.

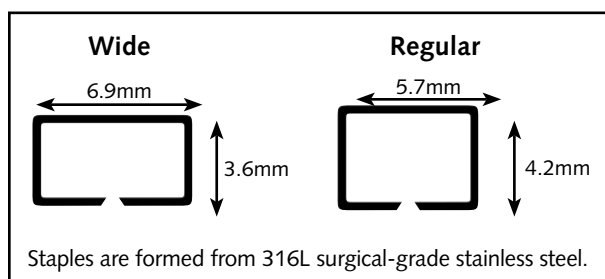

Reflex® One Skin Stapler

- Flexibility of movement.
- Fully Retracting Anvil permits staple release while moving forward or backward unlike competitive models which require you to "back off" before moving forward.
- Reflex® One releases the staple when the handle is fully depressed. Complete staple formation is confirmed by Audible and Tactile Feedback.

Reflex® TL Tissue Lift Skin Stapler

- Additional control.
- The Reflex® TL incorporates a Semi Retracting Anvil that only releases a staple when the handle is fully released.
- This Semi Retracting Anvil gives the TL its "Tissue Lift" which assists in approximating the incision for easier placement of the next staple.

Excellent Visibility

- Constant, visible indication of available staple count promotes more efficient use of supplies.
- Arrow guide for accurate placement.
- Open Nose Piece and low profile allow clear view of staple placement.

Superior Ergonomics

- Ergonomic design that helps prevent staple embedding and staves off the onset of Edema.
- Reflex®'s shape and balance fit any hand size for better feel and smoother action.

ORDERING INFORMATION

Product Description	Packaging	Reference #
Reflex® One Skin Stapler — Wide		
35 Wide	6/box	8535
25 Wide	6/box	8525
15 Wide	6/box	8515
Reflex® One Skin Stapler — Regular		
35 Regular	6/box	8035
25 Regular	6/box	8025
15 Regular	6/box	8015
Reflex® TL Skin Stapler with Tissue Lift — Regular		
35 Regular	6/box	8635
25 Regular	6/box	8625
15 Regular	6/box	8615
Reflex® TL Skin Stapler with Tissue Lift — Wide		
35 Wide	6/box	8735
25 Wide	6/box	8725
15 Wide	6/box	8715
Reflex® Extractor — Wide/Regular	12/box	3031
CONMED "Value" Extractor — Wide/Regular	12/box	3034

CONMED CORPORATION PRODUCT AREAS:

ARTHROSCOPY • ELECTROSURGERY • ENDOSCOPY • ENDOSURGERY • GASTROENTEROLOGY • INTEGRATED SYSTEMS • PATIENT CARE • POWERED INSTRUMENTS • PULMONOLOGY



525 French Road
 Utica, NY 13502
 (315) 797-8375
 Customer Service: 1-800-448-6506
 FAX: 1-800-438-3051
 International Sales: +1(315) 797-8375 • FAX: +1(315) 735-6235
 email: info@mail.conmed.com
 www.conmed.com

© 2009, CONMED Corporation, 04/2009, Control #12358 Rev. B



MBR 9010 Rev. 2