

AUTOCLAVE  
**PRIMA**



# Autoclave Prima

TECHNICAL EVOLUTION, GREAT DESIGN

## Our Galbiati Lab, 100% Italian

Department of “research and development”, created PRIMA, Class B **steam sterilizer** with an uninitiated design sense of aesthetic cleanliness combined with technology state-of-the-art and easy to use.

Solid, easy to use with an uncommon elegance, it's the new born of the Galbiati family, the first Italian production company of dental units.


## TECHNICAL SPECIFICATION

- **Rated voltage:** a.c.220V-230V, 50Hz
- **Rated power:** 18L 1500VA, 23L 1700VA
- **Material:** stainless steel (for medical)
- **Max. work pressure:** 2.5 bar Min. work pressure: -0.9 bar
- **Max. temperature:** 145 °C
- **Chamber volume:** 18L (Φ245×320mm) 23L (Φ245×450mm)
- **Loading size:** 18L (198×204×285) 23L (198×204×385)
- **Max. loading weight:** 18L (3.07kg/cm<sup>2</sup>)  
23L (3.21kg/cm<sup>2</sup>)
- **Working pressure/temperature:** 1.10~1.30bar /  
121°C~122°C; 2.10~2.30bar / 134°C~135°C
- **Water volume for one cycle:** 0.16L (min) 0.18L(max)

## PERFORMANCES

- Fast cycles (26min for HOLLOW 134°C)
- Limited water consumption
- Automatic water download
- 100% stainless steel
- Electronic water quality control
- Easy to position system
- GKE helix conformity test
- Integrated printer and usb data saving systems
- 24 months warranty

## DIMENSIONS

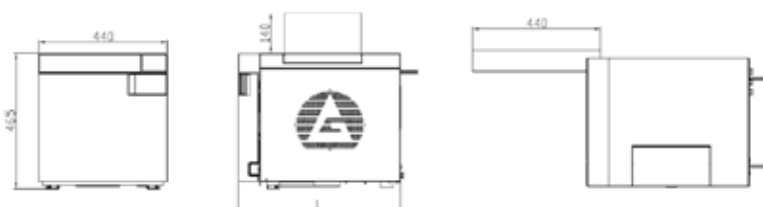

Item	18L Model	23L Model
Size of Closed Door	L615×465×440	VL715×465×440
Size of Open Door	L1005×605×440	L1105×605×440

## Button-control panel

Simple and intuitive

## Clean water tank

Stainless steel tank, indestructible

## Display

Cycle diagram display  
for live work control

## Chamber

4 storey tray chamber, usable  
also for surgical kits



## Quick connectors

Quick and easy water supply and  
drainage connections



## Main PCB board

Easy to reach PCB board for  
maintainance

# Galbiati

EXPERIENCE AND QUALITY SINCE 1935

“When a product does not work how it should, or our customers don’t like it, for us is a challenge. We feel that from our factory only good products, functional and accessible to all must be produced.

Our daily commitment is to satisfy our customer at 360° and be present before, during and after the purchase of Galbiati products, with the same passion and dedication”.



[www.galbiati.com](http://www.galbiati.com)  
[info@galbiati.com](mailto:info@galbiati.com)