



Argo Knotless[®]

Suture Anchor

Argo Knotless[®] Anchors feature Smart Tension Technology, ensuring consistent suture tension from insertion through deployment. Ready to eliminate unnecessary steps with a reproducible technique? Add in CONMED's Y-Knot[®] PRO RC All-Suture Anchors for a total cuff solution designed for simplicity.





BioBrace® Reinforced Bioinductive Implant
23 mm x 30 mm

The Complete Rotator Cuff Solution

When used in conjunction with Y-Knot® PRO RC All-Suture Anchors, Argo Knotless® Anchors provide surgeons with a complete, simple, and reliable solution for double-row rotator cuff repairs.

Treating a massive rotator cuff tear with poor tissue quality? Use Argo Knotless® Anchors in combination with CONMED's BioBrace® reinforced bioinductive implant for a strong, augmented repair.

- **Self-Punching**


A titanium tip configuration enables you to skip tunnel preparation and move straight to anchor insertion. Prefer to drill or punch? Choose our non-self-punching option.

- **Smart Tension Technology**

Set and maintain repair tension during anchor insertion, eliminating the risk of over-constraining the tissue.


- **Quick Release Driver**

After anchor insertion, the Argo Knotless® handle allows you to easily disengage the driver from the anchor, eliminating the need to backmallet or use other inefficient disengagement methods.

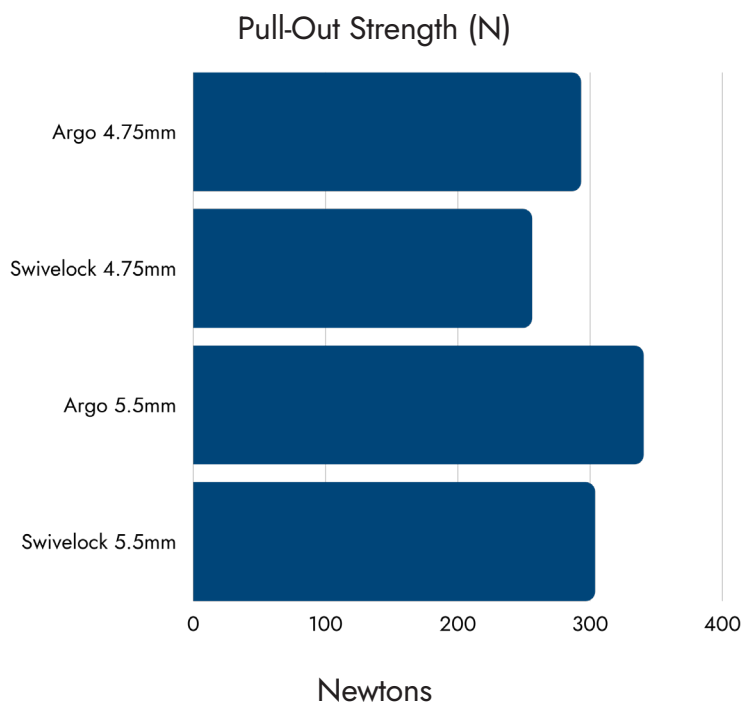

All anchors were loaded with 6 limbs of #2 Hi-Fi® sutures Istrom tester was used 40/12.5 pcf foam bone to simulate good bone quality.

Report code: PDD#1811996

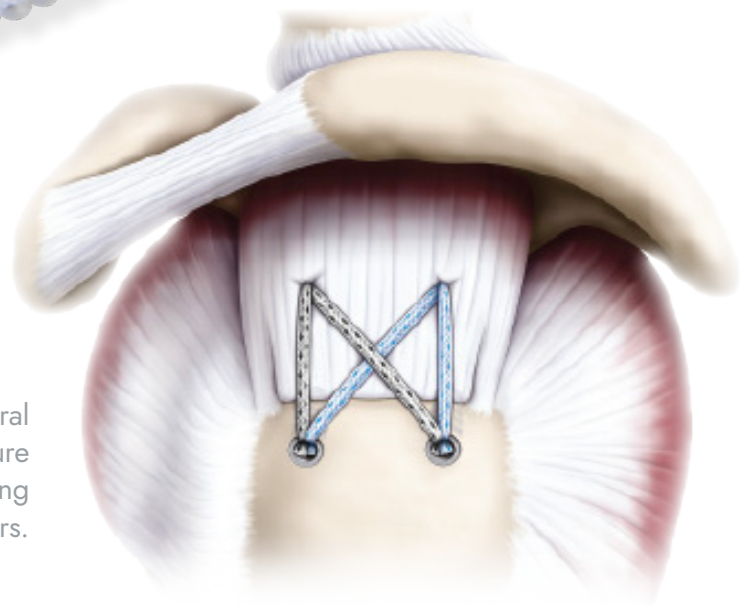

Argo Knotless® Anchors are the perfect lateral row counterpart to Y-Knot® PRO RC All-Suture Anchors – CONMED's most popular, best-selling medial row anchors.

Take the next step towards fewer steps in your technique and call **1.866.4.CONMED** to schedule an evaluation with your local CONMED Sales Representative, or visit www.CONMED.com/Argo for more information





Description

Catalog Numbers

Anchors

4.75mm Argo Knotless® Anchor with 1mm Hi-Fi® Ribbon	K475
5.5mm Argo Knotless® Anchor with 1mm Hi-Fi® Ribbon	K55
4.75mm Argo Knotless® Self-Punching Anchor with 1mm Hi-Fi® Ribbon	SPK475
5.5mm Argo Knotless® Self-Punching Anchor with 1mm Hi-Fi® Ribbon	SPK55

Instruments

4.75/5.5mm Argo Knotless® Disposable Drill Bit	K4755D
4.75/5.5mm Argo Knotless® Disposable Smooth Punch	K4755DSP
4.75/5.5mm Argo Knotless® Disposable Tapered Broaching Punch	K4755DTBP
4.75mm Argo Knotless® Disposable Tap	K475DTAP
5.5mm Argo Knotless® Disposable Tap	K55DTAP
4.75/5.5mm Argo Knotless® Reusable Smooth Punch	K4755RSP
4.75/5.5mm Argo Knotless® Reusable Tapered Broaching Punch	K4755RTBP
4.75mm Argo Knotless® Reusable Tap	K475RTAP
5.5mm Argo Knotless® Reusable Tap	K55RTAP

CONMED Corporation
11311 Concept Blvd
Largo, FL 33773

Toll Free: 1-866-4CONMED
International: 727-214-3000
www.CONMED.com

customerexperience@conmed.com
internationalorders@conmed.com