


RZ Medizintechnik GmbH
Unter Haslen 20 | D-78532 Tuttlingen
Phone: +49 (0)7462/9470-0
Fax: +49 (0)7462/9470-50
sales@rz-medizintechnik.com
www.rz-medizintechnik.com



MEMORY


Micro Instruments

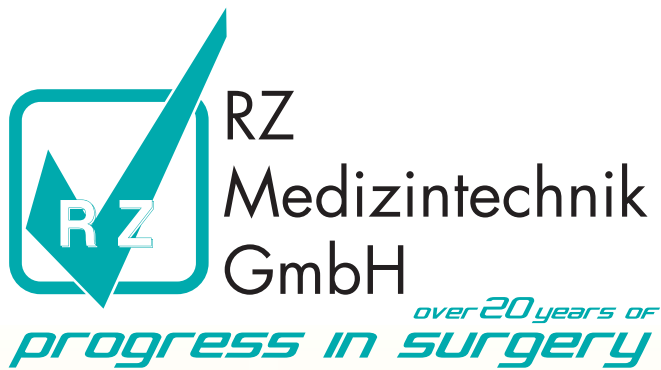

news
03.2014


MEMORY MICRO INSTRUMENTS

The RZ Memory Micro Instruments were developed for the use in sinuscopic and laryngoscopic operations. The **malleable shaft** and the **rotatable working ends** allow an easy adjustment of the instruments to meet any anatomical demands.

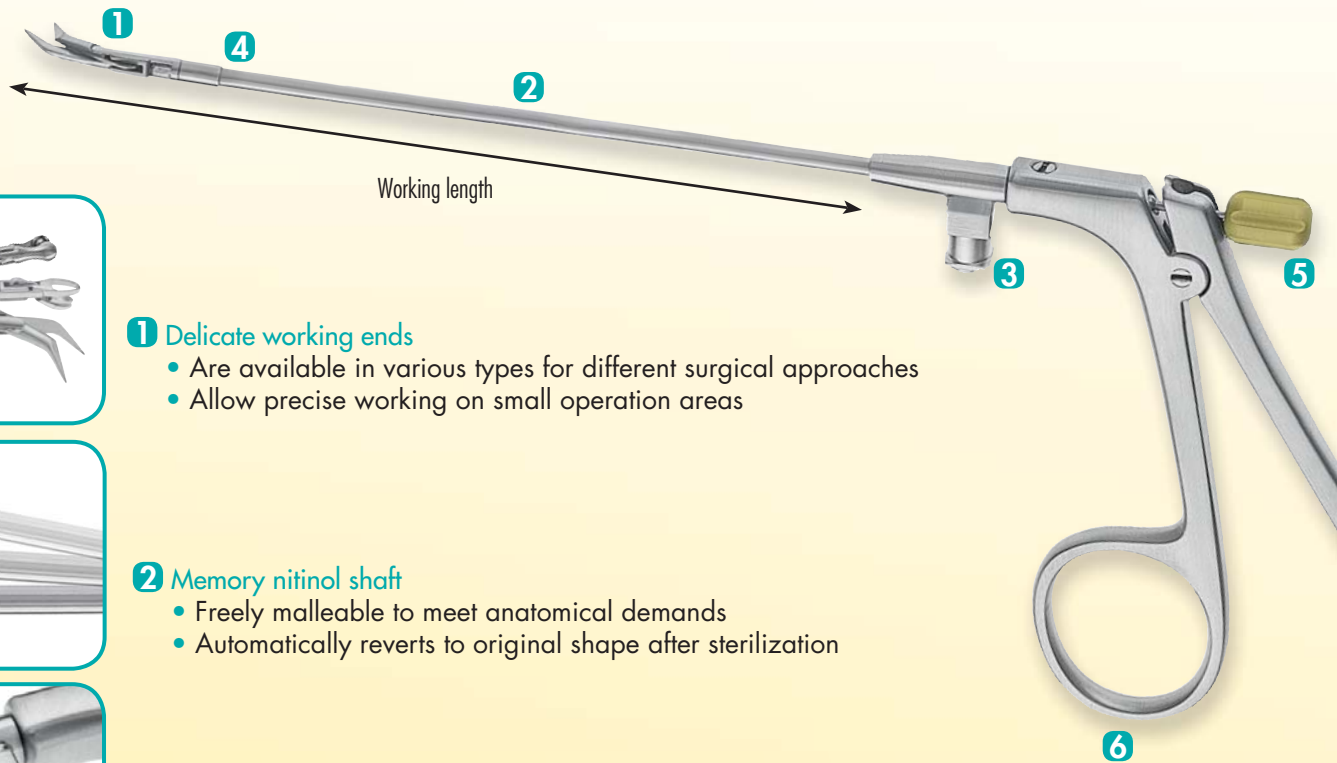

1 Delicate working ends

- Are available in various types for different surgical approaches
- Allow precise working on small operation areas

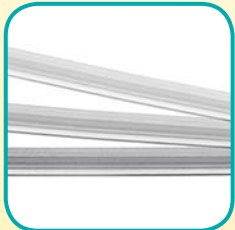

2 Memory nitinol shaft

- Freely malleable to meet anatomical demands
- Automatically reverts to original shape after sterilization



3 Irrigation connection

- Enables easy and efficient cleaning of internal components



4 360° rotatable working ends

- Ensure exact positioning by hand controlled rotary knob
- Allow unparalleled instrument positioning (in combination with malleable shaft)



5 Colour coated rotary knobs

- Ease identification of the instruments (working ends)



6 Lightweight ergonomic ring handle

- Reduces signs of fatigue and ensures comfortable working

Watch product video here! 



SCISSORS



Straight, sharp, single action

Sinuscopy
 | 252-800-000 WL 110mm
 | 252-800-140 WL 150mm
 Laryngoscopy
 | 100-800-230 WL 230mm



Straight, blunt, single action

Sinuscopy
 | 252-800-001 WL 110mm
 | 252-800-141 WL 150mm
 Laryngoscopy
 | 100-800-231 WL 230mm



Straight, sharp, double action

Sinuscopy
 | 252-800-002 WL 110mm
 | 252-800-142 WL 150mm
 Laryngoscopy
 | 100-800-232 WL 230mm



Curved, sharp, double action

Sinuscopy
 | 252-800-003 WL 110mm
 | 252-800-143 WL 150mm
 Laryngoscopy
 | 100-800-233 WL 230mm



Straight, sharp, double action

Sinuscopy
 | 252-800-004 WL 110mm
 | 252-800-144 WL 150mm
 Laryngoscopy
 | 100-800-234 WL 230mm



Curved, blunt, double action

Sinuscopy
 | 252-800-005 WL 110mm
 | 252-800-145 WL 150mm
 Laryngoscopy
 | 100-800-235 WL 230mm



45° angled, sharp, single action

Sinuscopy
 | 252-800-007 WL 110mm
 | 252-800-147 WL 150mm
 Laryngoscopy
 | 100-800-237 WL 230mm

GRASPING FORCEPS



Serrated, Ø 3mm

Sinuscopy

| 252-802-000 WL 110mm

| 252-802-140 WL 150mm

Laryngoscopy

| 100-802-230 WL 230mm

STRAIGHT FORCEPS



0.5mm tip

Sinuscopy

| 252-803-000 WL 110mm

| 252-803-140 WL 150mm

Laryngoscopy

| 100-803-230 WL 230mm



0.9mm tip

Sinuscopy

| 252-803-001 WL 110mm

| 252-803-141 WL 150mm

Laryngoscopy

| 100-803-231 WL 230mm

CURVED FORCEPS



Curved, fine

Sinuscopy

| 252-801-000 WL 110mm

| 252-801-140 WL 150mm

Laryngoscopy

| 100-801-230 WL 230mm



Curved, wide

Sinuscopy

| 252-801-001 WL 110mm

| 252-801-141 WL 150mm

Laryngoscopy

| 100-801-231 WL 230mm

BLAKESLEY FORCEPS



Sharp, double action

Sinuscopy
| 252-805-000 WL 110mm
| 252-805-140 WL 150mm
Laryngoscopy
| 100-805-230 WL 230mm



Through cutting, single action

Sinuscopy
| 252-807-000 WL 110mm
| 252-807-140 WL 150mm
Laryngoscopy
| 100-807-230 WL 230mm



Sharp, Ø 2mm

Sinuscopy
| 252-804-002 WL 110mm
| 252-804-102 WL 150mm
Laryngoscopy
| 100-804-232 WL 230mm



Sharp, Ø 3mm

Sinuscopy
| 252-804-003 WL 110mm
| 252-804-103 WL 150mm
Laryngoscopy
| 100-804-233 WL 230mm



Sharp, Ø 5mm

Sinuscopy
| 252-804-005 WL 110mm
| 252-804-105 WL 150mm
Laryngoscopy
| 100-804-235 WL 230mm

SPOONS



© 10/2015 Copyright by RZ Medizintechnik GmbH - Unter Haslen 20 · 78532 Tuttlingen · Germany. All rights reserved.
Modifications of technical specification or design are subject to change without prior notice.