

Unter Haßlen 20 · D-78532 Tuttlingen
Phone: +49 (0)7462/9470-0
Fax: +49 (0)7462/9470-50
e-mail: sales@rz-medizintechnik.com
<http://www.rz-medizintechnik.com>



OSCILLATING SAW autopsy & plaster



14/08

Complete range of saw blades for Autopsy and Plaster saws



200-115-200 Autopsy saw, 230V, detachable, autoclavable
 200-115-201 Autopsy saw, 110V, detachable, autoclavable
 For use in pathology departments (post-mortem)
 Saw head can be removed for steam sterilization

- Slim and lightweight design
- Silent motor and gear
- Automatic load-compensation over the complete oscillating field
- Safety ON/OFF switch
- Solid and maintenance-free design
- Radio/Audio-interference-suppressed

Complete Set includes Autopsy Saw, round saw blade 65mm Ø,
 round saw blade 50mm Ø,
 pair of wrenches and
 hexagonal key.



200-115-100 Plaster saw, 230V
 200-115-101 Plaster saw, 110V

200-118-100 Vacuum Cleaner 230V
 200-118-101 Vacuum Cleaner 110V

200-115-500 Carry case with foam inlay



Round saw blade 65mm Ø




Round saw blade 50mm Ø



Round saw blade 45mm Ø




Post Mortem saw blade 65mm Ø

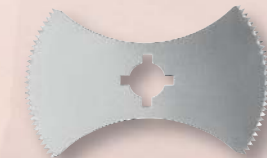

Segmented saw blade 65mm Ø



Depth saw blade 8mm wide

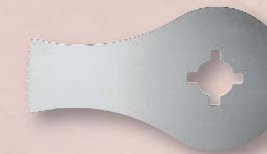

Depth saw blade 20mm wide



200-115-005 Eccentric Disc to limit depth of cut



200-115-006 Pair of wrenches



200-115-007 Hexagonal key

Plaster / Autopsy

Synthetic casts

PTFE coated
(non stick)

200-115-165

200-116-165

200-117-165

200-115-150

200-116-150

200-117-150

200-115-145

200-116-145

200-115-265

200-116-265

200-115-365

200-116-365

200-115-108


200-116-108

200-115-120

200-116-120

TECHNICAL DETAILS OF ALL RZ AUTOSPY AND PLASTER SAWS

- 500 Watt power (Autopsy Saw), 210 Watt power (Plaster Saw)
- 6500 to 24000 oscillations, infinitely variable (Autopsy saw),
10000 to 19000 oscillations, infinitely variable (Plaster Saw)
- Fully automatic power-compensation above the whole oscillations field
- Slim motor edition with silent motor 52 – 77 dB(A) depending on to
the selected oscillations
- Latest electronic system with various protection and security functions
- Special positioning of the motor, protected by plaster dust
- Robust construction for highest loading capacity by continuous operation
- Safety ON / OFF switch
- Weight only 1,25 kg without cable (1,5 kg with cable) for Autopsy Saw
Weight only 1,00 kg without cable (1,2 kg with cable) for Plaster Saw
- Vibrations- and maintenance-free
- Regulation from outside of the manual operating area
- 5 meter cable with EU-connection, 230 Volt version (Autopsy Saw)
3 meter cable with EU-connection, 230 Volt version (Plaster Saw)
3 meter cable with US-connection, 110 Volt version
- Double protection isolation according to class II, medical use of typ BF
- Radio/Audio-interference-supressed
- CE conformity
- EN 60601-1, EN 60601-1-2, IEC 60601-1, IEC 60601-1-2 inspected

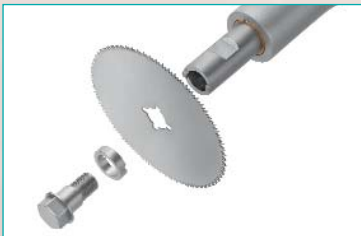

Universal docking system for all saw blades



Easy and quick change of saw blades



6500 – 24000 oscillations,
infinitely variable (Autopsy Saw)
10000 – 19000 oscillations,
infinitely variable (Plaster Saw)

