


Quality time for better care

"Quality time for better care" is the brand promise of Terumo Medical Care Solutions

HOSPITAL CARE

Terufusion™

IV Administration Set-
PT+PU400WY
Product Information

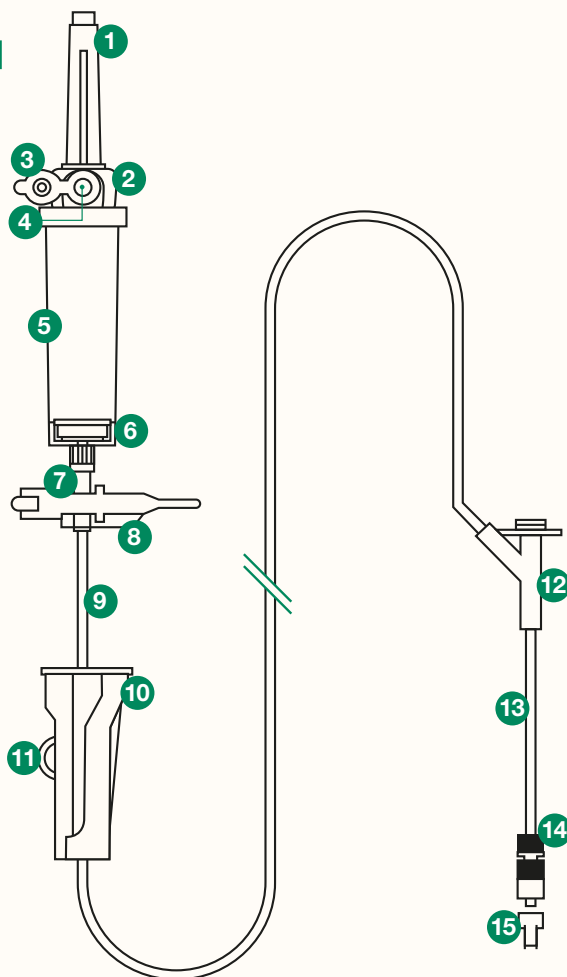

Description		Raw Materials
1	Protective cap	Polyethylene
2	Spike	ABS
3	Air filter cap	PVC
4	Air filter	Glass fiber
5	Drip chamber	PVC
6	Filter	ABS and Nylon
7	Cover tube	PVC
8	Anti-Free Flow clamp	Polyoxymethylene
9	Tubing (Ø 3.15 x 2000 mm)	PVC
10	Flow regulator housing	ABS
11	Flow regulator wheel	ABS
12	Injection point (Y-Type)	ABS + Polyisoprene
13	Tubing (Ø 3.15 x 200 mm)	PVC
14	Rotating Luer-Lock	ABS
15	Luer cap with filter	LDPE

Does not contain components made of DEHP or natural rubber latex.



Sterilization

With ethylene oxide to a Sterility Assurance Level of at least 10⁻⁶.



Storage

- Shelf life of 5 years
- Batch number and expiry date are printed on the label on the pouch
- Do not store between 5°C and 35°C

IV administration set
for use with Terumo Terufusion™
infusion pumps series¹



A clear view on infusion with a comprehensive portfolio of solutions



Packaging

Single Unit Pouch

Unit Box 26 pieces

Dimensions 275 x 230 x 135 mm, Gross Weight 1.28 kg

Shipping Carton 250 pieces

Dimensions 683 x 430 x 342 mm, Gross Weight 13.90 kg



Precautions

- Do not detach the Anti Free Flow Clip (AFF Clip) from IV administration set.
- We recommend to switch the tubing section used in the pump after maximum 24 hr to another position as the flow rate may deviate from the set value due to deformation of the tubing section in the pump after long period of use.
- Dispose of safely after use to avoid risk of infection.



Conformity

EN ISO 13485, EN ISO 11135, EN ISO 15223-1, EN ISO 8536-4, EN ISO 8536-8, EN ISO 10993-1, ISO 14644-1, EN ISO 11607-1, EN 868-5

Highly accurate
infusion combination of
Terumo pump & sets¹



Product Details

Product code	PT+PU400WY
Spike	2-way spike
Vented	Yes
Drip chamber	20 drops/mL
Liquid filter	15 µm
Roller clamp	Orange with white wheel
Anti-Free Flow (AFF) clamp	18.4 ml
Y-site	With needle port
Priming volume (empty drip chamber)	18.1 ml
Total length	2380 mm



Does not contain components made of DEHP or natural rubber latex.



Regulatory Status

CE Certification – Class IIa medical device

Notified Body

TÜV SÜD 0123

Manufacturer

Poly Medicure Ltd. (India)

EC Representative

OBELIS SA, Boulevard General Wahis 53, 1030 Bussels, BELGIUM

EC Distributor

TERUMO EUROPE N.V.

Interleuvenlaan 40, 3001 Leuven, Belgium

Your Terumo partner



TERUMO EUROPE N.V.

Terumo Corporation Terumo Europe N.V.
+81 3 6742 8500 +32 16 38 12 11

www.terumo-europe.com

Reference

1 See data on file

All brand names are trademarks or registered trademarks owned by TERUMO CORPORATION, its affiliates, or unrelated third parties. Published by Terumo Europe N.V.